



An easy spectrophotometric acid-base titration protocol for dissolved organic matter

Marawit Tesfa, Aline Dia, Malou Ollivier, Iván David Osorio Leon, Remi
Marsac

► To cite this version:

Marawit Tesfa, Aline Dia, Malou Ollivier, Iván David Osorio Leon, Remi Marsac. An easy spectrophotometric acid-base titration protocol for dissolved organic matter. *MethodsX*, 2022, 9, pp.101721. 10.1016/j.mex.2022.101721 . insu-03664468

HAL Id: insu-03664468

<https://insu.hal.science/insu-03664468>

Submitted on 11 May 2022

HAL is a multi-disciplinary open access archive for the deposit and dissemination of scientific research documents, whether they are published or not. The documents may come from teaching and research institutions in France or abroad, or from public or private research centers.

L'archive ouverte pluridisciplinaire **HAL**, est destinée au dépôt et à la diffusion de documents scientifiques de niveau recherche, publiés ou non, émanant des établissements d'enseignement et de recherche français ou étrangers, des laboratoires publics ou privés.



Distributed under a Creative Commons Attribution 4.0 International License

An easy spectrophotometric acid-base titration protocol for dissolved organic matter

Marawit TESFA , Aline DIA , Malou OLLIVIER ,
Iván David OSORIO LEON , Rémi MARSAC

PII: S2215-0161(22)00102-9
DOI: <https://doi.org/10.1016/j.mex.2022.101721>
Reference: MEX 101721

To appear in: *MethodsX*

Received date: 29 September 2021
Accepted date: 28 April 2022

Please cite this article as: Marawit TESFA , Aline DIA , Malou OLLIVIER ,
Iván David OSORIO LEON , Rémi MARSAC , An easy spectrophotometric acid-base titration protocol for dissolved organic matter, *MethodsX* (2022), doi: <https://doi.org/10.1016/j.mex.2022.101721>



This is a PDF file of an article that has undergone enhancements after acceptance, such as the addition of a cover page and metadata, and formatting for readability, but it is not yet the definitive version of record. This version will undergo additional copyediting, typesetting and review before it is published in its final form, but we are providing this version to give early visibility of the article. Please note that, during the production process, errors may be discovered which could affect the content, and all legal disclaimers that apply to the journal pertain.

© 2022 Published by Elsevier B.V.

This is an open access article under the CC BY license (<http://creativecommons.org/licenses/by/4.0/>)

An easy spectrophotometric acid-base titration protocol for dissolved organic matter

Marawit TESFA^{1*}, Aline DIA¹, Malou OLLIVIER¹, Iván David OSORIO LEON¹, Rémi MARSAC¹

¹University of Rennes, CNRS, Géosciences Rennes - UMR 6118, Rennes, France

*Corresponding author: marawit.tesfa@univ-rennes1.fr


Abstract

UV-vis spectrophotometric acid-base titration can characterize dissolved organic matter (DOM) acid-base properties. However, it requires incremental pH adjustment, which make the procedure time consuming and results are subjected to dilution effect. This study brings forth a new approach, referred as the “buffer method” for pH adjustments by using carefully selective pH-buffers to adjust the pH. This, statistically validated method minimizes the pH adjustment time and lightens the laboratory work load. Chemical product cost associated with this novel method is slightly increased as compared to the previous approach, due to the necessity to use pH-buffers.

- Buffer method: Acid-base titration by using buffer for pH adjustment
- Buffer method validated by statistical means
- Rapid, reliable and economical method

Keywords: acid-base titration; UV-vis spectrophotometry; dissolved organic matter

Graphical abstract

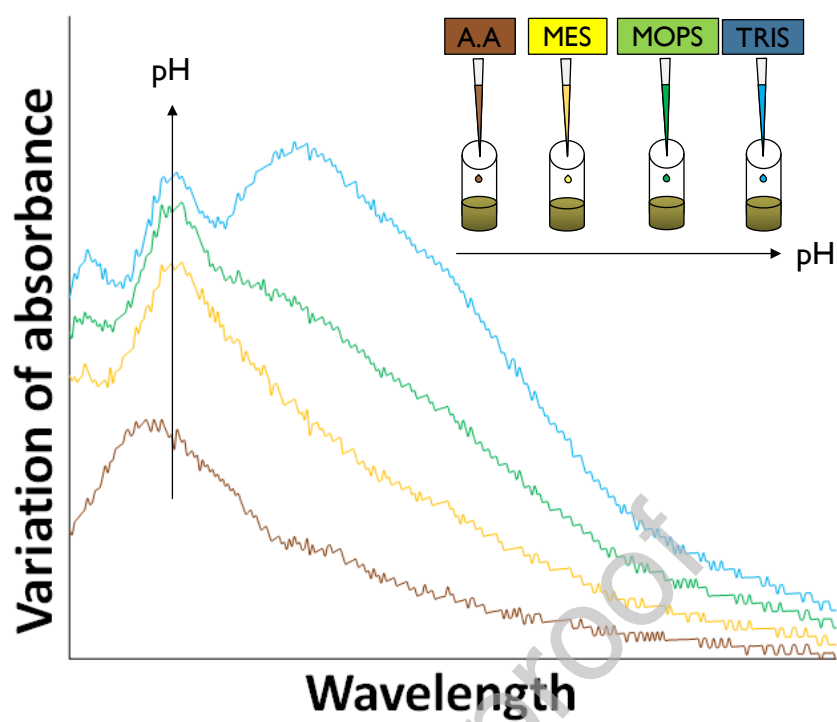

<u>Specifications Table</u>	
Subject area	Environmental sciences
More specific subject area	Wet chemistry
Method name	“Buffer method”
Name and reference of original method	Original method has no name in the literature, it was named in this article “Incremental method” References of original method Dryer, D.J., Korshin, G.V., Fabbicino, M., 2008. In Situ Examination of the Protonation Behavior of Fulvic Acids Using Differential Absorbance Spectroscopy. Environ. Sci. Technol. 42, 6644–6649. https://doi.org/10.1021/es800741u Gao, Y., Yan, M., Korshin, G.V., 2015. Effects of Ionic Strength on the Chromophores of Dissolved Organic Matter. Environ. Sci. Technol. 49, 5905–5912. https://doi.org/10.1021/acs.est.5b00601 Janot, N., Reiller, P.E., Korshin, G.V., Benedetti, M.F., 2010. Using Spectrophotometric Titrations To Characterize Humic Acid Reactivity at Environmental Concentrations. Environ. Sci. Technol. 44, 6782–6788. https://doi.org/10.1021/es1012142 Liu, S., Benedetti, M.F., Han, W., Korshin, G.V., 2020. Comparison of the properties of standard soil and aquatic fulvic and humic acids based on the data of differential absorbance and fluorescence spectroscopy. Chemosphere 261, 128189. https://doi.org/10.1016/j.chemosphere.2020.128189
Resource availability	All resources used in this study can be acquired from Sigma Aldrich and the International Humic Substances Society

Dissolved organic matter (DOM) is ubiquitous in aquatic systems and plays a key role in the environment as it is at the crossroad of multiple biogeochemical cycles. In particular, because of its capacity to form complexes with metal contaminants, the quantification of chemical entities that drive

DOM-metal ion interaction is highly required (Benedetti et al., 1995; Marsac et al., 2021, 2014; Milne et al., 2003, 2001; Tipping, 1998). These are mainly acid-base functional groups, which can be determined by acid-base titrations. Amongst various approaches, UV-vis spectrophotometric acid-base titration allows the analysis of DOM acid-base properties at environmentally relevant concentrations (Dryer et al., 2008; Gao et al., 2015; Janot et al., 2010; Liu et al., 2020). This can be used, for instance, to monitor changes of DOM deprotonation–protonation properties in water treatment processes (Yan et al., 2016).

The principle of spectrophotometric acid-base titration is to progressively increase the pH of the DOM solution with consecutive additions of small volumes of base. The absorbance of the solution is measured by spectrophotometry after each step of pH increment. Despite the simplicity of the method, it has not yet been automated, making the laboratory workload meticulous and time-consuming. In addition, as a consequence of all of these consecutive steps, this experimental procedure presents a significant dilution effect. This commonly used method will be referred to as the “incremental method”.

To lighten the laboratory workload and minimize the dilution effect with the incremental method, we developed an alternative to the pH stabilization. Based on studies that have used buffers to keep the pH of DOM solutions constant without affecting the optical properties of the DOM (Gao et al., 2015; Liu et al., 2020), we proposed to generalize the buffer use to a wider range of pH, necessary for acid-base titration. The newly developed “buffer method” referred to as such further in this study, consists in dividing a DOM mother solution into aliquots, adjusting these aliquots pH individually with carefully selected buffers and measuring aliquots absorbance independently. Furthermore, this new protocol can be also automated using an autosampler connected to the UV-vis spectrophotometer. The “buffer method” presents an automatable, faster experimental setup yielding results with minimum dilution effects for UV-vis spectroscopic acid-base titration of DOM.

Reagents. Suwannee River Fulvic Acid (SRFA, ref n° 3S101F) was acquired from the International Humic Substance Society (IHSS). This fulvic acid sample was selected because it is most commonly used and therefore, widely characterized. Potentiometric titration conducted by Ritchie and Perdue, (2003) indicate that SRFA sample contains 12.00 mg L^{-1} of carboxylic groups and 1.48 mg L^{-1} of phenolic groups. The pH-buffer salts (Table 1), background salt (NaCl), acid (HCl) and base (NaOH) solutions were purchased from Sigma Aldrich. The pH-buffers are some of the widely used “Good’s buffers” and acetic acid/acetate, which are available at relatively low cost (about 1€/g or below).

Instruments. Dissolved Organic Carbon (DOC) concentrations were determined by Shimadzu TOC-L carbon analyzer. The pH was monitored using HANNA, HI5221 pH meter. The pH meter was calibrated with standard HANNA buffers (4.01, 7, 10.1) and, new calibrations were performed for each buffer solution preparation. Absorbance was recorded, in a 1-cm-wide quartz cuvette, between 200 and 800 nm, with a UV-vis Shimadzu UV2600 spectrophotometer.

Solution preparation. A stock of DOM mother solution (1 g L^{-1}) was prepared by dissolving 250mg of organic matter powder into a total 250 mL volume of milliQ water and base (0.1 M NaOH). Complete dissolution was guaranteed by maintaining solution pH at 9 with NaOH for 24 hours with adjustments after 3, 6 and 12 hours. The mother solution was diluted to 25 mgC L^{-1} and the ionic strength (IS) to 0.01 M with background electrolytes (NaCl, 1 mol L^{-1}). Buffer mother solutions were made by dissolving buffer salts in 40 mL ultrapure water. The pH of the buffer solutions was adjusted to the pKa value of the buffer or within a range of $\text{pKa} \pm 1$ by using acid (HCl, 1 mol L^{-1}) or base (NaOH, 1 mol L^{-1}), then verified and readjusted if necessary after 3, 6, 12 and 24 hours. The final volume was raised to 50mL with ultrapure water, leading to a final buffer concentration of 1 mol L^{-1} .

Incremental titration protocol. The “incremental method” (Janot et al., 2010) is illustrated in Figure 1a and was carried out as follows: a 50 mL volume of the DOM mother solution was placed in a 100 mL plastic polypropylene beaker with a magnetic stirrer and under N_2 bubbling throughout the titration process. The pH of the DOM solution in the beaker was adjusted to 3 using HCl (0.1 mol L^{-1}). The quartz cuvette was filled with 3 mL of the solution and the UV-vis spectrum was recorded. The 3 mL were dropped back to the beaker to avoid any variation of the volume of the titrated DOM solution. The pH of the titrated solution was slightly raised to 3.2 by manual addition of a small volume of NaOH (0.1 mol L^{-1}). After pH stabilization and recording, the quartz cuvette was rinsed with the DOM solution then analyzed by spectrophotometry, and then put back into the beaker. Spectra were measured every ca. 0.2 pH units between 3 and 9. The volume of each base addition was

recorded to later be considered for dilution calculations. Considering the pH stabilization and the pH modification steps (0.2 unit), this experimental process takes at least three to four hours for a single sample, most of the time being dedicated to pH stabilization process.

Titration protocol using pH buffers. The principle of the newly developed “buffer method” is to split the DOM solution into several aliquots (13 in the present study) and to adjust their pH using pH buffers. The pH of two aliquots was adjusted by HCl addition at pH 3 and 3.4. The pH of the 11 other aliquots was fixed using the buffer solutions, details of the buffer are listed in Table 1. The aliquots contained a buffer solution at a concentration of $10^{-3} \text{ mol L}^{-1}$. Absorbance was measured once the pH of each aliquot stabilized. This protocol was illustrated in Figure 1b. Considering 10 mins for pH adjustment per sample (with HCl and buffer), 15 extra mins for pH stabilization, and manual absorbance measurements, this experimental setup only takes at most one to one and a half-hour. The UV-vis spectra of $10^{-3} \text{ mol L}^{-1}$ buffer solutions in 10 mM NaCl in the absence of DOM were also recorded.

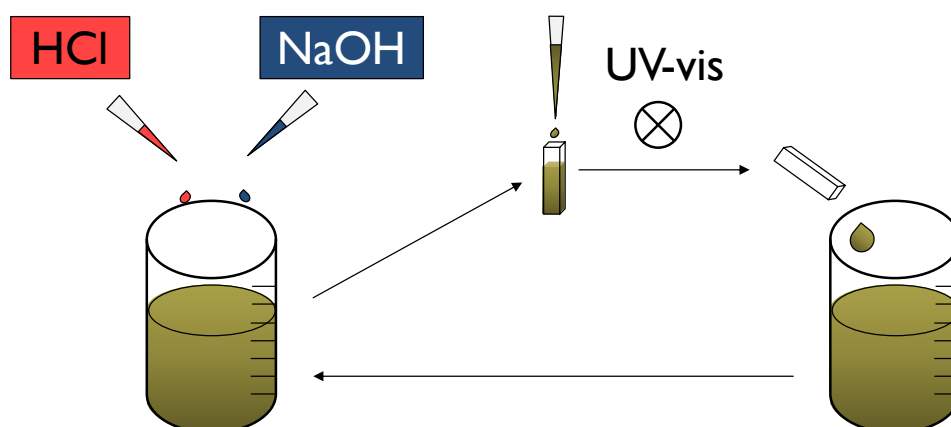
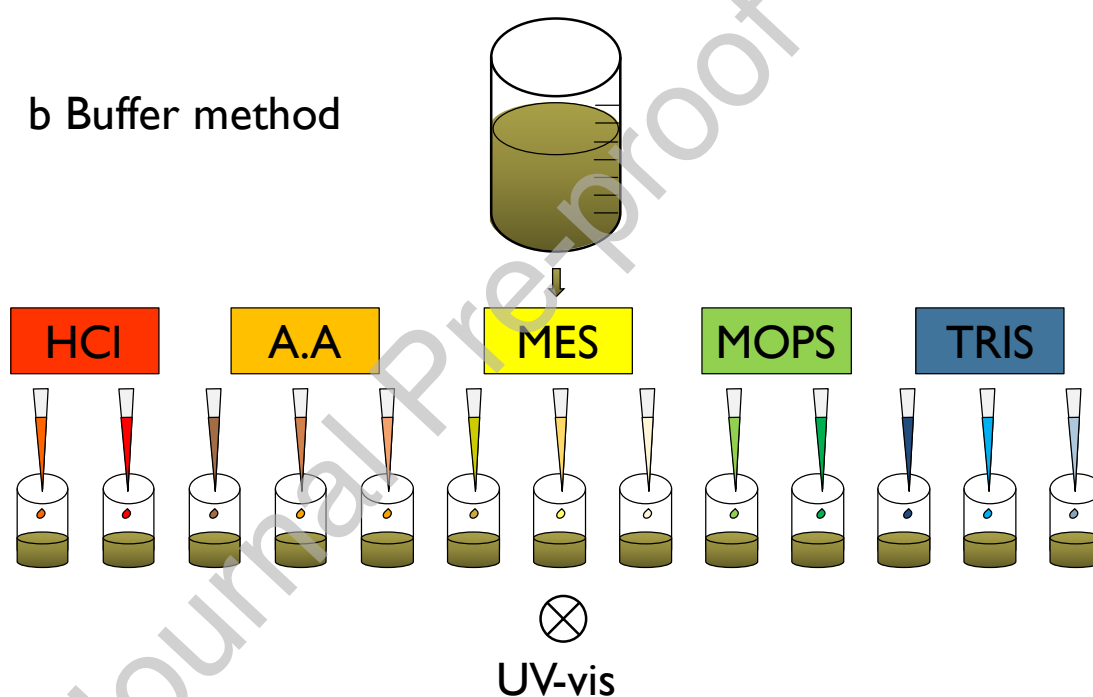
a Incremental method**b Buffer method**

Figure 1. Scheme of the experimental setup for the “incremental method” (a) and the “buffer method” (b).

Table 1: Buffer solution properties used in this study

Buffer	Chemical formula	pKa	pH
TRIS	$C_4H_{11}NO_3$	8,1	9.06
			8.62
			8.09
MOPS	$C_7H_{14}NaO_4S$	7,15	7.44
			7.18
			6.44
MES	$C_6H_{13}NO_4S$	6,15	6.11
			5.51
			4.91
AA	$C_2H_4O_2$	4,1	4.48
			4.00

Application example

Absorbance of the buffer solutions. The optical property of each buffer was investigated in a blank buffer solution in 10 mM NaCl, i.e. in the same conditions as the titrated DOM solution. Results are shown in Figure 2. For wavelengths lower than 240 nm, a high absorbance was recorded for all of the buffer solutions, which is also observed for the background electrolyte solution (NaCl only). For some buffers, a slight absorption peak is observed within the 240 and 290 nm range with a maximum absorbance of $0.014 \text{ cm}^{-1}\text{g}^{-1}\text{L}^{-1}$, which is very small and represents less than 1% (0.7% on average) of the absorbance of the corresponding samples containing 25 mgC L^{-1} DOM (Figure 2). As the main wavelengths of interest, for DOM characterization, are centered around 270 and 370 (Janot et al., 2010), it is safe to consider that the buffer molecules do not contribute significantly to the absorbance of the DOM-containing samples.

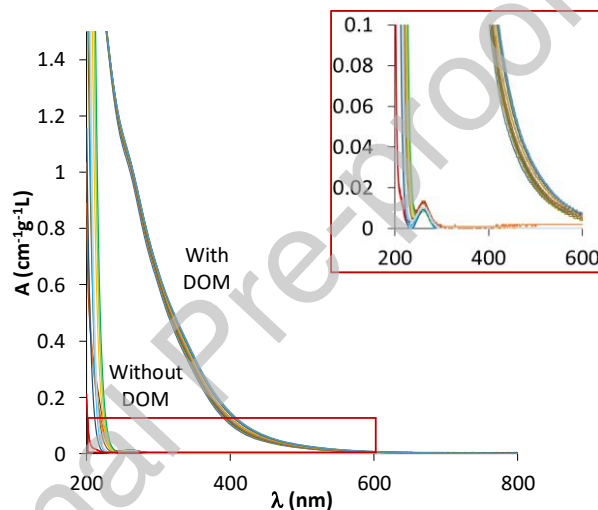


Figure 2: Absorbance spectra of buffer solutions (1 mM) in 10 mM NaCl, in the absence (colored lines) or the presence of 25 mgC L^{-1} DOM (black lines).

Spectrophotometric titration results. The absorbance of DOM samples slightly varies with pH (Figure 2) (Janot et al., 2010). For visual ease, the absorbance spectra (A) were presented as the variation of absorbances (ΔA) with respect to a reference spectrum as a function of the wavelength (λ) for both incremental and buffer methods in Figure 3a and 3b, respectively. Their corresponding raw absorbances are presented in inserts in each figure. The sample at $\text{pH} = 3$ was taken as a reference to determine ΔA values:

$$\Delta A_{\lambda} = A_{\lambda,i} - A_{\lambda,\text{pH}=3} \quad (1)$$

The absorbance measurement, with both incremental and buffer methods, presents a typical DOM absorbance spectra; the absorbance monotonously and progressively decreases with wavelength and increases with pH (Janot et al., 2010). With both methods (Figures 3 a and b) two peaks appear. The first one, centered around 270, is major at all pH values, whereas the second one, centered around 370, is smaller and appears as a shoulder of the first one only for circumneutral to alkaline solutions (i.e for $\text{pH} > 7$). Comparing the two sets of ΔA , none difference stands out suggesting that the selected buffers do not affect the optical acid-base properties of the tested sample within the considered wavelength range.

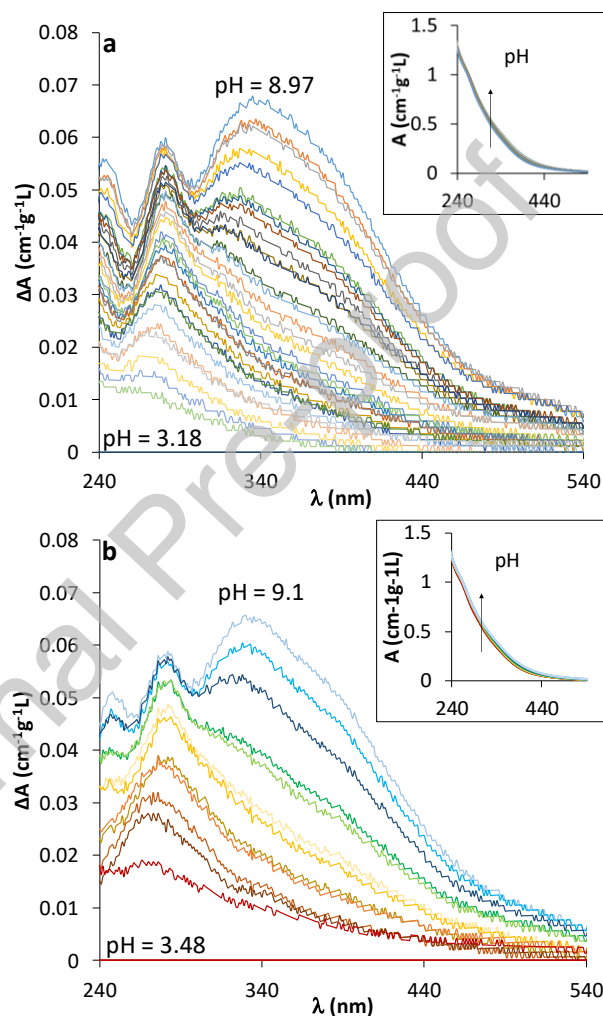


Figure 3: Relative absorbance (ΔA) measurement with the incremental (a) and buffer (b) methods. Inserts in each figure are the corresponding raw absorbance (A).

Validation. To verify that the use of buffer does not significantly impact the measured wavelength, a statistical validation was performed. The datasets of the incremental and the buffer methods were compared, at each wavelength, through mean and regression comparisons. The means of the incremental and the buffer ΔA were statistically compared by Student's t-test ($p < .001$) for each wavelength. After hypothesis verification, none of the wavelengths present a significant difference between the data drawn from the incremental and the buffer methods. Since comparing the means for

each wavelength set could omit other information such as the variation, two linear models were developed to explain the correlation between absorbance and pH: one for the “incremental method” and other for the “buffer method”. The slope of the two models were then compared by a Student’s t-test. None significant difference was underlined for all of the considered wavelengths ($p > .001$) when comparing the two models representing the sample pH stabilization method. This result suggests that the sample preparation has no significant impact on the absorbance variation (Figure 4).

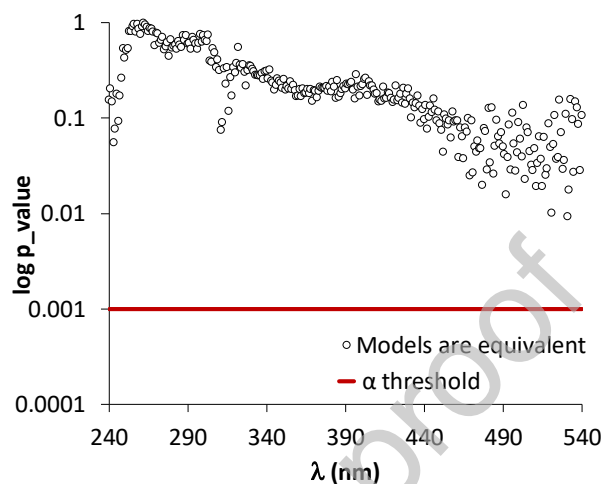


Figure 4: Evolution of p-values versus the wavelength in the regression analysis between the incremental and buffer methods. P-values correspond to a slope comparison using a Student’s t-test. The significance threshold α is 0.001.

Final remark This study brings forth a new experimental setup to characterize the acid-base properties of DOM by spectrophotometry acid base titration. Using the “buffer method”, an automatable, faster and reliable protocol, characterizing DOM at environmentally relevant concentrations becomes less tedious. Considering a coupling between a spectrophotometer and an automatic sampler, this method would allow an automation of the acid-base titration by UV-vis spectrophotometry, thus reducing the laboratory work period to a minimum. The disadvantages of this new method are only minor as it requires (i) additional chemicals (i.e. pH-buffers) associated with limited cost and (ii) the preparation of pH-buffer stock solutions, but that can be used for numerous analyses. Multiplying the characterization of unknown DOM will help create and feed a data bank that will allow a better understanding of the acid-base properties of DOM and, indirectly, of its reactivity towards metal contaminants in the environment.

Declaration of Competing Interest

The authors of this article declare that they have no conflict of interests.

Acknowledgments

This work was supported by the C-FACTOR project funded by ANR (project number ANR-18-CE01-0008). Through the support of the GeOHeLiS analytical platform of Rennes University, this publication is also supported by the European Union through the European Regional Development Fund (FEDER), the French ministry of Higher Education and Research, the French Region of Brittany and Rennes Metropole.

References

- Benedetti, M.F., Milne, C.J., Kinniburgh, D.G., Van Riemsdijk, W.H., Koopal, L.K., 1995. Metal Ion Binding to Humic Substances: Application of the Non-Ideal Competitive Adsorption Model. *Environ. Sci. Technol.* 29, 446–457. <https://doi.org/10.1021/es00002a022>
- Dryer, D.J., Korshin, G.V., Fabbicino, M., 2008. In Situ Examination of the Protonation Behavior of Fulvic Acids Using Differential Absorbance Spectroscopy. *Environ. Sci. Technol.* 42, 6644–6649. <https://doi.org/10.1021/es800741u>
- Gao, Y., Yan, M., Korshin, G.V., 2015. Effects of Ionic Strength on the Chromophores of Dissolved Organic Matter. *Environ. Sci. Technol.* 49, 5905–5912. <https://doi.org/10.1021/acs.est.5b00601>
- Janot, N., Reiller, P.E., Korshin, G.V., Benedetti, M.F., 2010. Using Spectrophotometric Titrations To Characterize Humic Acid Reactivity at Environmental Concentrations. *Environ. Sci. Technol.* 44, 6782–6788. <https://doi.org/10.1021/es1012142>
- Liu, S., Benedetti, M.F., Han, W., Korshin, G.V., 2020. Comparison of the properties of standard soil and aquatic fulvic and humic acids based on the data of differential absorbance and fluorescence spectroscopy. *Chemosphere* 261, 128189. <https://doi.org/10.1016/j.chemosphere.2020.128189>
- Marsac, R., Banik, N.L., Marquardt, C.M., Kratz, J.V., 2014. Stabilization of polynuclear plutonium(IV) species by humic acid. *Geochim. Cosmochim. Acta* 131, 290–300. <https://doi.org/10.1016/j.gca.2014.01.039>
- Marsac, R., Catrouillet, C., Davranche, M., Bouhnik-Le Coz, M., Briant, N., Janot, N., Otero-Fariña, A., Groenenberg, J.E., Pédrot, M., Dia, A., 2021. Modeling rare earth elements binding to humic acids with model VII. *Chem. Geol.* 567, 120099. <https://doi.org/10.1016/j.chemgeo.2021.120099>
- Milne, C.J., Kinniburgh, D.G., Tipping, E., 2001. Generic NICA-Donnan Model Parameters for Proton Binding by Humic Substances. *Environ. Sci. Technol.* 35, 2049–2059. <https://doi.org/10.1021/es000123j>
- Milne, C.J., Kinniburgh, D.G., van Riemsdijk, W.H., Tipping, E., 2003. Generic NICA-Donnan Model Parameters for Metal-Ion Binding by Humic Substances. *Environ. Sci. Technol.* 37, 958–971. <https://doi.org/10.1021/es0258879>
- Ritchie, J.D., Perdue, E.M., 2003. Proton-binding study of standard and reference fulvic acids, humic acids, and natural organic matter. *Geochim. Cosmochim. Acta* 67, 85–96. [https://doi.org/10.1016/S0016-7037\(02\)01044-X](https://doi.org/10.1016/S0016-7037(02)01044-X)
- Tipping, E., 1998. Humic Ion-Binding Model VI: An Improved Description of the Interactions of Protons and Metal Ions with Humic Substances. *Aquat. Geochem.* 4, 3–47. <https://doi.org/10.1023/A:1009627214459>
- Yan, M., Dryer, D., Korshin, G.V., 2016. Spectroscopic characterization of changes of DOM deprotonation–protonation properties in water treatment processes. *Chemosphere* 148, 426–435. <https://doi.org/10.1016/j.chemosphere.2016.01.055>

Journal Pre-proof